

Technical Product Report

For Research Use Only; Not for use in Diagnostic Procedures

Product Description: Seraseq® ctDNA Reference Material WT

Material No: 0710-0720 Batch Number: 10387904

Material Description: Human genomic DNA from the reference cell line, GM24385.

Date of Manufacture: 27-Dec-18 Expiration Date: 27-Dec-20

Concentration (Qubit dsDNA BR Assay): Nominal value: 20 ng/mL; Average measured value after extraction using Qiagen QIAamp Circulating Nucleic Acid Kit: 24.6 ng/mL


Volume: 5 mL

Storage: 2-8 °C

Digital PCR testing using BioRad QX200™ Droplet Digital™ PCR System:

Gene ID	COSMIC Identifier	Amino Acid Change	Average AF%
EGFR	COSM6240	p.T790M	0.01
EGFR	COSM6224	p.L858R	0.00
EGFR	COSM6225	E746_A750 delELREA	0.00
EGFR	COSM6252	p.G719S	0.00

Approval:

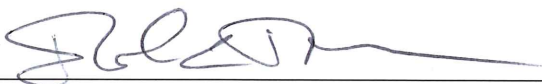

Prepared By

10-Jan-19

Date