

AccuPlex™ SARS-CoV-2 Reference Material

ONGOING PERFORMANCE MONITORING FOR ROUTINE PATIENT TESTING

HIGHLIGHTS

COMPATIBLE WITH CDC AND WHO CONSENSUS SEQUENCES FOR SARS-COV-2

REAL MAMMALIAN VIRUS ENABLES FULL PROCESS QC FOR RT-PCR BASED VIRAL DISEASE ASSAYS

NON-INFECTIOUS AND REPLICATION DEFICIENT; ENSURES SAFETY DURING ASSAY DEVELOPMENT AND END-USER QC APPLICATIONS

AccuPlex™ SARS-CoV-2 Reference Material Kit is designed to measure day-to-day performance of the assay. The product contains positive materials directed against the published CDC and WHO consensus sequences (Table 1), and negative materials targeting the human RNase P gene. The materials are formulated for use with *in vitro* diagnostic test methods that detect SARS-CoV-2 virus, the causative agent of COVID-19 disease.

Product Description	Material Number	Sequence Coverage	Pack Size	Concentration
AccuPlex™ SARS-CoV-2 Reference Material Kit	0505-0126	ORF1a, RdRp, E, N	Positive 5 x 1.5 mL	5,000 copies/mL
			Negative 5 x 1.5 mL	5,000 copies/mL (RNase P)


Not for In Vitro Diagnostic Use. Research Use Only.

SARS-COV-2 ACCUPLEX SOLUTION:

- Non-infectious and replication deficient, enables safe handling in contrast to viral samples
- Fully-extractable with a real viral protein coat; superior to “naked” transcribed RNA
- Long shelf life: 2 year stability at 2 - 8°C

Table 1: Regions included in AccuPlex Coronavirus (SARS-CoV-2) Positive Reference Material

Region	Included Nucleotides
ORF1a	417..1899
	3094..3360
RdRp	13291..13560
	14700..15950
	18577..19051
E (Envelope)	25801..28200
N (Nucleocapsid)	27952..29873

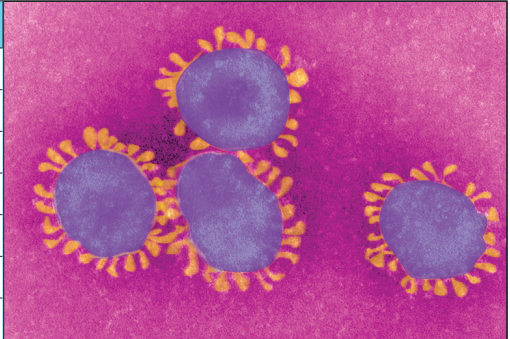

Numbers refer to sequence of GenBank NC_045512.2

To learn more about AccuPlex SARS-CoV-2 Quality Solutions visit: www.seracare.com/SARSCoV2 or contact customer service at 800.676.1881

ABOUT LGC SERACARE

TRUSTED SUPPLIER
TO THE DIAGNOSTIC
TESTING INDUSTRY
FOR OVER 30 YEARS

HIGH-QUALITY
CONTROL PRODUCTS,
RAW BIOLOGICAL
MATERIALS, AND
IMMUNOASSAY
REAGENTS

INNOVATIVE TOOLS
AND TECHNOLOGIES
TO PROVIDE
ASSURANCE IN
DIAGNOSTIC ASSAY
PERFORMANCE AND
TEST RESULTS

FOR MORE
INFORMATION, PLEASE
VISIT OUR WEBSITE:
WWW.SERACARE.COM

LGC SeraCare's AccuPlex recombinant material serves as a true full process molecular control for your diagnostic assays. Compatible with multiplexed RT-PCR and NGS-based assays, AccuPlex custom recombinant materials are constructed with a replication-deficient mammalian virus, producing a safe, non-infectious material (Figure 1). With a protein coat and lipid bilayer, these mammalian virus-based reference materials resemble the complexity of virus targets found in true patient samples.

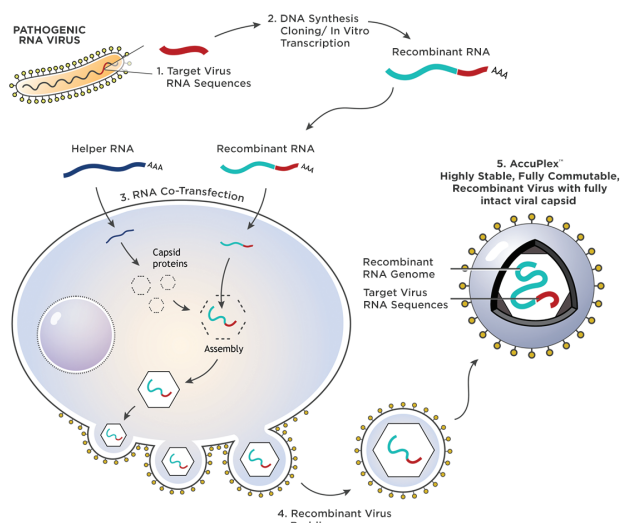


FIGURE 1: 1) RNA sequence from the pathogenic virus of interest is chosen. 2) DNA synthesis and cloning occur to produce the recombinant RNA. 3) Recombinant RNA and helper RNA are co-transfected into the mammalian cells, allowing the encapsulation of recombined RNA. 4) Exocytosis of the mature enveloped non-infectious and replication deficient RNA virus containing the assay target RNA sequence of interest

MOST 'PATIENT SAMPLE-LIKE' MATERIAL EVALUATES ENTIRE WORKFLOW

Unlike RNA-spiked buffer, or technologies that package viral RNA into a bacteriophage, the AccuPlex recombinant closely resembles the wild-type mammalian pathogenic virus. This enables the release of the viral genome at a similar rate to the wild-type virus during the nucleic acid sample preparation process. AccuPlex recombinant material mimics a real patient sample in your workflow, serving as a full-process control for your assay.

ACCUPLEX SOLVES ASSAY DEVELOPMENT CHALLENGES

If you're developing diagnostics for emerging viral diseases and have the challenge of including safe, non-infectious controls in your test kit, partner with LGC SeraCare's talented R&D team to produce your custom AccuPlex recombinant virus material (DNA or RNA-based). Utilizing your sequences of interest and product specifications, we will develop a custom solution which meets your unique requirements.

ORDERING INFORMATION

To place an order, or learn more about our SARS-CoV-2 Quality Solutions, please contact us at +1.508.244.6400 and 800.676.1881 or email customerservice@seracare.com



AccuPlex™ is a trademark of SeraCare Life Sciences, Inc.
© 2021 SeraCare Life Sciences, Inc. All rights reserved.
MKT-00610-01